

### SYSTEM SPECIFICATIONS

Laser Medium	Nd:YAG Crystal
Wavelength	1064nm / 532nm
Operating Mode	1064nm / 532nm / Andy Toning [PTP] / Multi Toning / Long [FR]
Beam Profile	JTS Top Mode
Pulse Energy	Max. 1,600mJ [ 1064nm / Andy Toning / Multi Toning ] Max. 500mJ [ 532nm ] Max. 2,200mJ [ Long(FR) ]
Pulse Rate	1 - 15 Hz
Cooling System	Closed cycle water to air heat exchanger
Pulse Width	5ns - 10ns [ 1064nm / 532nm / Andy Toning / Multi Toning ] 300us [ Long(FR) ]
Spot Size	Zoom Handpiece : 2 - 10mm
Power Input	230V-, 50/60Hz, 3KVA
Dimensions	650 X 850 X 270mm [ W X H X D ]
Weight	65kg

# Q-Switch Nd YAG



**ANDY PRO**  
ND:YAG LASER








A dramatic experience where skin turns white like snow blossom

# SNOW BLOSSOM

\* LASER TONING

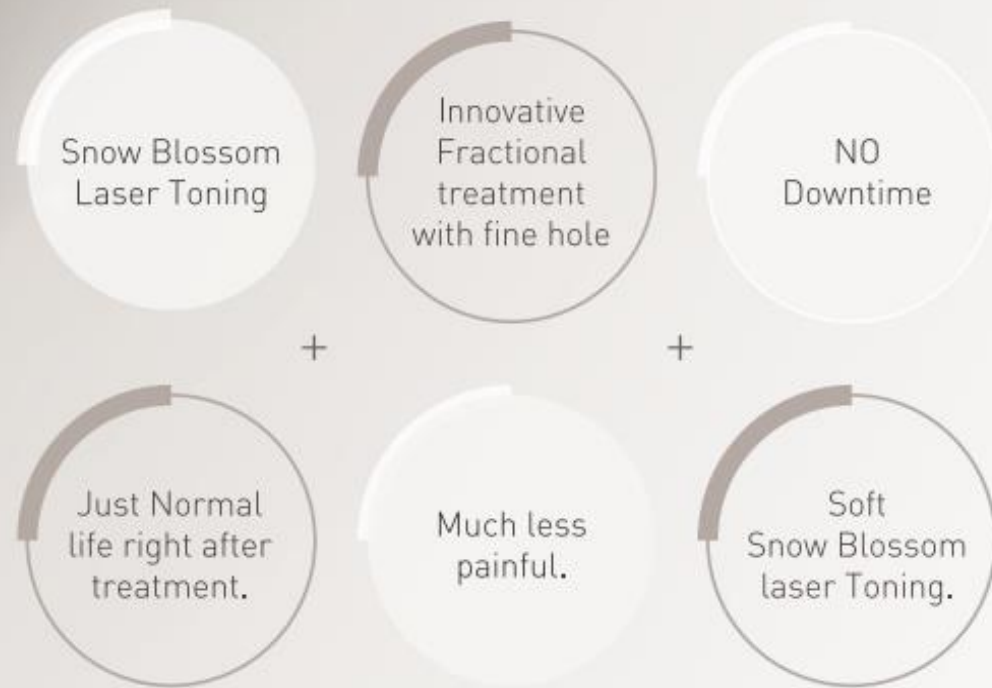

My skin turns white like a white snow blossom

# SNOW BLOSSOM

\* LASER TONING

ANDY Pro is providing the unique laser toning, Snow Blossom.

Conventional laser toning can't do but Snow Blossom could deliver the best protocol for the new generation of laser toning with the innovative Fractional treatment.

- Treatment Laser Toning
- No Anesthesia
- The period of treatment within 10 mins.
- No Downtime.

**ANDY PRO**  
ND:YAG LASER



## Indications

- > Feel the innovative Fractional laser Toning
- > The new generation of Tattoo Removal
- > PIH and mild treatment of pigmentation
- > To get the better skin texture and tone.

Freckle & blemish	Seborrheic Keratosis	Caf'e-Au-Lait spot
Tattoo (Red)	ABNOM/Hor's Nevus	Tattoo (Black, blue, brown)
Eye Line	Eye Brows	An-Dy Toning
Inflammatory Acne	Erythema	Enlarged Pores



## JTS HANDPIECE



Zoom handpiece [2-10mm]

2,3,4,5,6,7,8,9,10mm size hand piece is available in zoom hand piece. The sensor recognize the spot size and the full 10.2 inch monitor shows the spot size and alarm if the spot size is set as different size.



Collimated handpiece [7mm]

7mm size of hand piece can help treatment of toning easily and accurately. With special optics inside of hand piece, energy loss is decreased and more focus laser energy on target area even though the laser hand piece is not contacted on the face tightly.



1064 M.L.A Fractional H.P

Spot size is 7mm and >121 Holes. Various spot sizes and shapes are available with fractional handpiece and enables large area's treatment.


1064 D.O.E Fractional H.P

Beam size is 5X5mm and 49 Holes. Various spot sizes and shapes are available with fractional handpiece and enables large area's treatment.

532 D.O.E Fractional H.P

Beam size is 4X4mm and 81 Holes. Various spot sizes and shapes are available with fractional handpiece and enables large area's treatment.


## 1064 & 532 Pulse mode of ANDY Pro

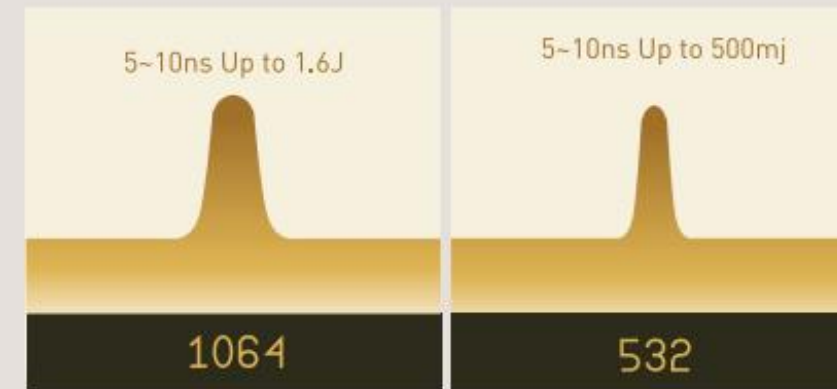

With short duration, in this mode Andy can support high power in a moment, It breaks melanic pigment from deep dermis to epidermis thorouly


532nm wavelength is shorter than 1064nm and its penetration depth area is superficial epidermis area. Q532 mode treat epidermal pigment lesions effectively

## Fractional Beam Mode

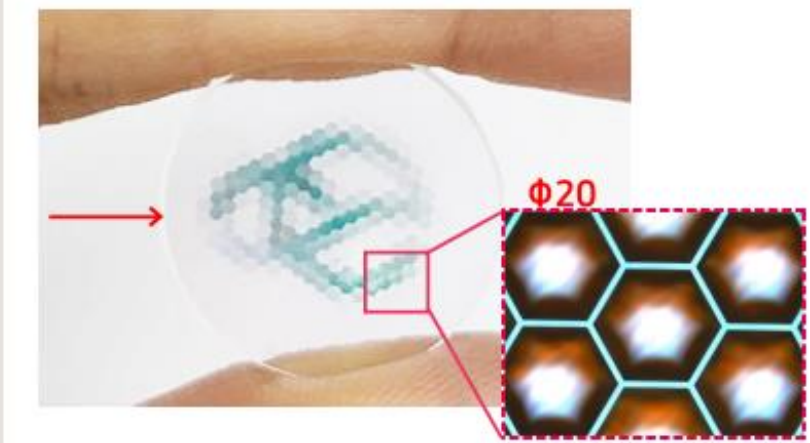

Combination effect of toning treatment by fractional hand piece


## Andy Toning / Long Peel mode

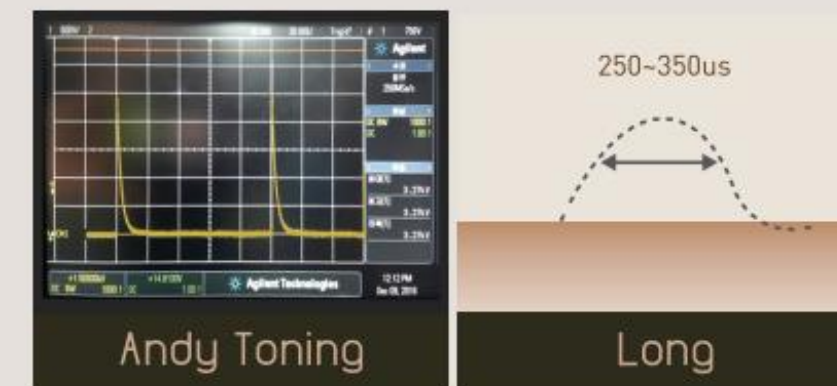

After first pulse is radiated on the skin, twin pulse is radiated within 100us of short interval and it increases effdecive of photo thermal effect and breaks melasma safely

With short duration, in this mode Andy can support high power in a moment, It breaks melanic pigment from deep dermis to epidermis thorouly

## JTS Heating 532nm Laser

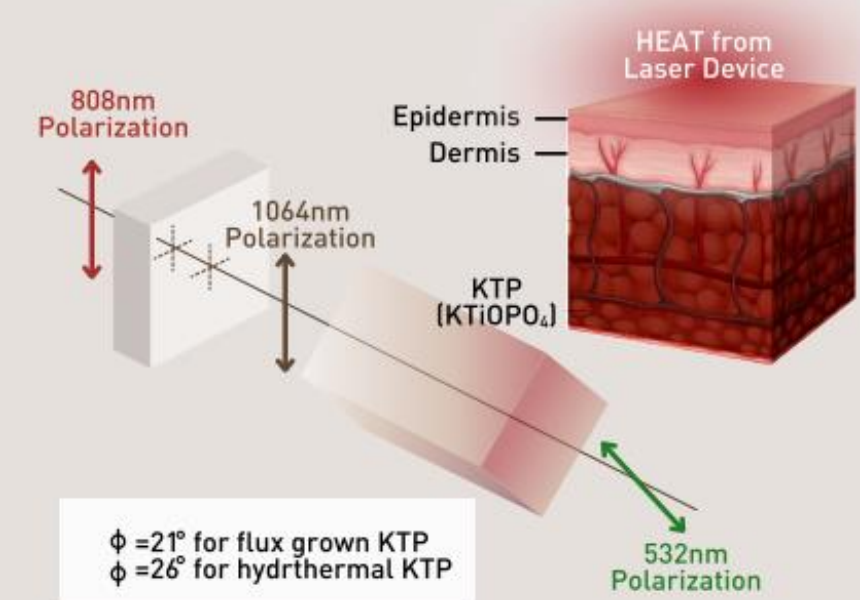

## JTS Beam Technology

